


Massager Factory Store

Dropshipping Contact WhatsApp:
+86 189 2273 5243

USER MANUAL

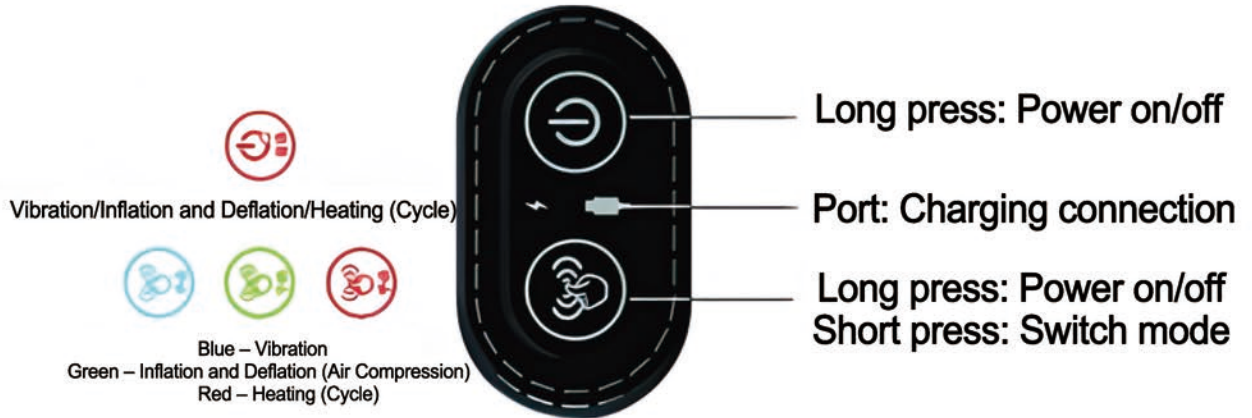

Detailed Instructions

1. Power Connection

Keep the device connected to power during use (use a 5V/2A adapter if necessary).

Once it is powered, the device performs an automatic self-test (the green light remains steadily on).

Dropshipping: Contact us via WhatsApp: +86 189 2273 5243

2. Starting the Device

To turn the device on or off, press and hold the button for 2 seconds.

3. Function Modes (Manual mode – requires button presses)

- Blue light: Vibration massage
- Green light: Airbag compression
- Red light: Heating function

4. RECOMMENDED AUTOMATIC MODE

Continuous heat + automatic switching between vibration and air compression.
Hands-free operation.

Manual Mode Option

Short button press to switch between functions (only one active at a time).

5. Safety Instructions

Use for a maximum of 30 minutes per session to avoid overheating.
Unplug after use.