


90X120MM 骑马钉 单色

SMART HAND MASSAGER

User Manual

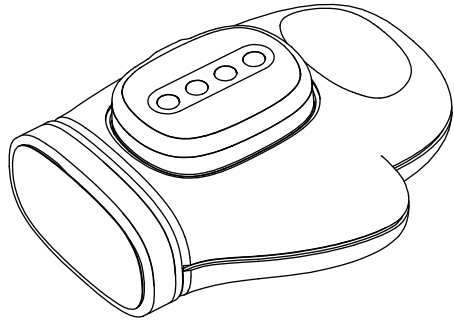

Please read the manual carefully before using the product and keep it properly.
The product appearance is subject to the actual product.

Precautions

When using this product, please read all instructions in this manual carefully, especially the precautions. Violating the instructions may damage the machine or even cause danger.

Note

1. Do not use for purposes other than hand massage.

2. Do not use with other therapeutic devices at the same time.

3. Do not let young children use or play.

4. Do not use when sleeping.

5. People who cannot express their wishes or cannot move independently should not use it.

6. This product is for one-handed use, do not use it with both hands at the same time.

7. When using, do not wear hard objects such as jewelry or nail decorations around the massage area.

8. Do not use it in places with heavy humidity such as bathrooms, and do not let liquids such as water come into contact with the product.

9. The device has a heating surface, and people who are not sensitive to heat must pay attention when using it.

10. The default working time of this product is 15 minutes per time. If you need to use it again, you need to turn it on again. It is recommended that the single massage time should not exceed 30 minutes to avoid excessive massage.

Prohibited people

1. People who have had bone correction or surgery on their hands are prohibited from using the machine.

2. People with severe peripheral nerve disorders and sensory impairments.

3. People with skin injuries.

4. People with acute pain such as osteoporosis, sprains, muscle strains, etc.

5. People who cannot operate the machine by themselves.

Usage and safety instructions

1. Please choose according to the strength you can bear and make sure it is close to the skin.

2. It is recommended to use the charging cable and standard adapter provided by the machine for charging.

3. Do not clean the machine with water.

4. When using it for the first time, or the battery power is depleted after the product has been shelved for a long time, please follow the instruction manual to charge it to ensure that it can work properly.

5. Do not use it with wet hands.

Precautions

Maintenance and storage

1. Cleaning

Wipe with a soft cloth, do not wet the main unit, do not wipe the main unit with an overly wet towel to prevent water from seeping into the main unit, and do not use thinners, gasoline or other volatile liquids to wipe.

2. Storage and storage

Place the instrument out of reach of children, away from moisture, high temperature, direct sunlight and water splashes. After use, it is recommended to store the product in the original box. If it is not used for a long time, it is recommended to charge it once every 30 days or so to avoid the battery triggering over-discharge protection, which may cause abnormal charging or failure to charge.

3. Maintenance

Do not disassemble, repair or replace parts by yourself, which may cause machine failure. If there is a problem with the product, please turn it off in time, do not continue to use it and contact the company's maintenance service center.


Product Parameters

Charging voltage	5V
power	5W
Battery capacity	2000mAh


Instructions

Product Introduction

Function buttons

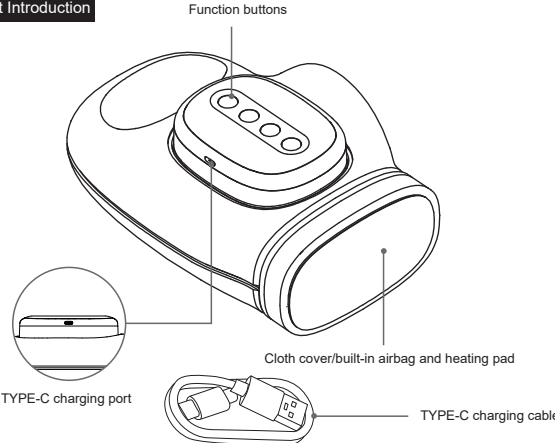

Use an adapter that complies with the standard, 5V=1A/5V=2A (Prepare Yourself)

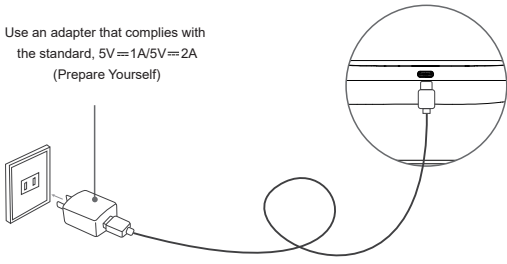


Preparation before use

Please check whether the product has power. If the power is too low and the product automatically shuts down or cannot be turned on, please charge it before using it.

1. Make sure your hands are clean and do not get them wet. Do not wear jewelry.

2. Put your fingers into the correct holes in order.

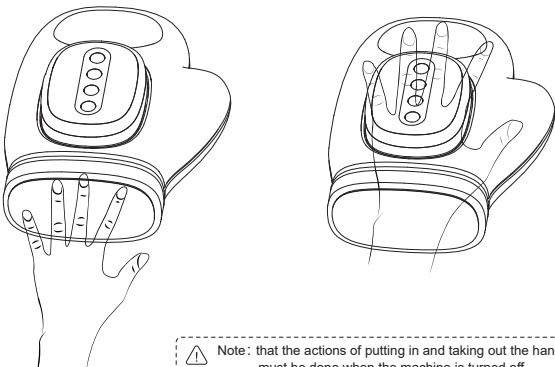

Note: that the actions of putting in and taking out the hand must be done when the machine is turned off.

Instructions

Introduction to left and right hand use

This product is for one-handed use, please do not use it with both hands at the same time.

Left hand

Right hand

Product Introduction

Heat button

Mode key

Velocity Keys

Power button

1. Power button

Long press "Power button" to turn on/off. When the battery is fully charged, the green light will be on; low battery: the red light will flash; charging: the red light will flash slowly; fully charged: the green light will be on

2. Hot compress button

Short press "Heat button" to switch between low/medium/high/on/off hot compress. The default setting is no heating. Low-grade hot compress, green light on; medium-grade hot compress, yellow light on; high-grade hot compress, red light on; off hot compress, light off.

3. Mode button

The default mode is one when powered on. Short press "Mode key" to switch the mode (mode 1-mode 2-mode 3). Mode 1, green light on; Mode 2, blue light on; Mode 3, red light on.

4. Velocity Keys

Short press "Velocity Keys" to switch between 1-5 strength levels, 1 strength levels are the default when powered on. Level 1, green light on; Level 2, blue light on; Level 3, yellow light on; Level 4, purple light on; Level 5, red light on.

Function	Air pressure	Heat	Time
Mode 1	Finger-palm-wrist (comprehensive)	/	15min
Mode 2	Finger-palm-wrist (ladder)	/	15min
Mode 3	Finger-palm-wrist (alternating)	/	15min

Battery Disassembly


Battery disassembly

Step 1: Open the host decoration

Step 2: Remove the screws

Step 3: Take out the host

Step 4: Remove the battery

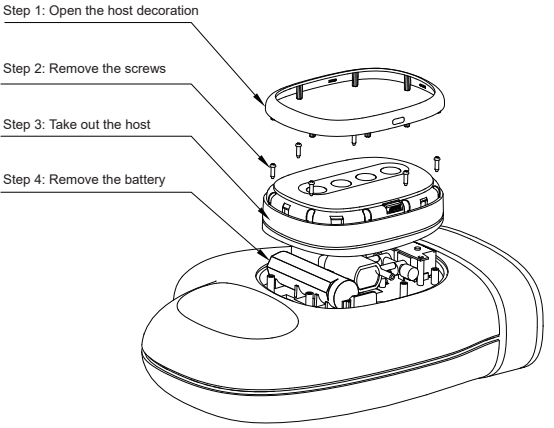

The main unit (including lithium battery) and accessories of this product have reached the end of their service life. Please do not discard them at will, otherwise electronic waste will be generated and pollute the environment. Before discarding the product, the battery must be removed from the product, before removing the battery, the product must be powered off, the battery should be disposed of safely when removed.

Troubleshooting

Troubleshooting

If you encounter problems when using the product, please refer to the following tips and solutions. If the product is still in a faulty state or cannot be used, please turn off the power. Users should not disassemble or repair the product by themselves. Disassembly or repair by yourself and damage to the product will cause the product warranty to be invalid, and the company will not be responsible for any consequences caused by this situation.

problem	Possible causes/solutions
Unable to start the product	Please confirm whether the battery is low, please confirm that the On-Off button has been pressed.
Cannot charge or no charging light	Please check whether the TYPE-C charging cable is properly connected to the charging port of this device;
The product suddenly stops while working.	Check if the product is low on power; The time set by the automatic timing device has passed; If the problem persists, do not continue to use it and contact the company's repair service center;

Danger

The battery used in this product is a special rechargeable battery and cannot be used in other products. Do not charge the rechargeable battery that has been removed; Do not throw it into fire to cause excessive heat; Do not pierce it with nails, impact it with gravity, or disassemble it for modification; Do not directly connect the positive and negative poles of the battery with metal; Do not store it with necklaces, hairpins, etc. Do not charge, use, or place it near a stove or in the hot sun. The above situations may cause heating, fire, rupture, etc.

Warning

Please place the removed rechargeable battery out of reach of infants and young children to prevent accidental swallowing; if accidental swallowing occurs, please seek medical attention immediately;
When the battery fluid leaks, do not touch the liquid with bare hands, and handle it as follows:
When the battery fluid enters the eyes, there is a risk of blindness. Do not rub the eyes directly with your hands. Immediately wash with clean water and seek medical attention;
When the battery fluid gets on the body or clothing, it may cause skin inflammation or injury. Please rinse with clean water. If you still feel uncomfortable, seek medical attention immediately;
The product can only be disassembled when it is discarded. Failure to do so may result in fire, electric shock, injury, etc.