

IPL laser hair remover User Manual




User guide

IPL Laser hair remover is used to remove excess hair and permanently reduces hair growth. Do it 1-3 times a week, after 4 weeks, can effectively remove excess hair. Suitable for body parts: face (only for below mouse), underarms, arms, legs, bikini lines, etc.

Principle

IPL laser hair remover are personal care appliances, for safety reasons, please read this user manual carefully before first time usage.

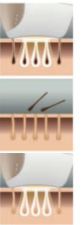

Notice

Read and follow all safety warnings and instructions carefully before use.


IPL laser hair removal will flash at the treatment site, do not look directly at the flash.

Tanned or exposed sunlight skin before, please do not use IPL laser hair remover.

Avoid skin exposure in the sun for 4 weeks before using IPL laser hair remover.

Do not use IPL laser hair remover near your eyes.

If your skin is naturally dark, do not use IPL laser hair removal.

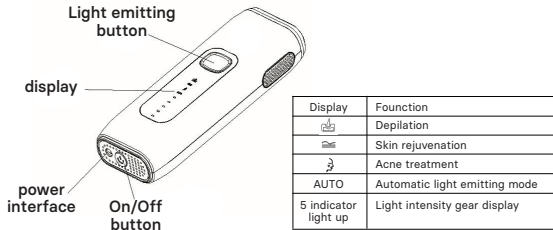

Parameter of product	
Product number: AI06	Power consumption: 36W
Input voltage: 24V	Current consumption: 1.5A

Contains in package

1. IPL laser hair remover machine
2. AC adapter
3. User manual
4. Protective glasses
5. Scrape knife

preparing for use

Before use, thoroughly clean the skin, shave the hair of the clean hair removal area (less than 1mm), ensure that the skin is free of hair, completely dry, no oily substance residue (perfume, lotion, body fluids, sunscreen and other residues) to ensure that light energy is guided to the root in the most effective Note: When using with IPL laser, there is no shaving or unclean skin, and various adverse conditions may occur during lighting. .
01. The dust in the hair, which may absorb light energy, may cause pain or redness of the skin, fork hair or dust contaminates the light window, preventing normal light may lead to reduced product life.
02. The hairs on the skin may burn and severely burn the skin on the surface.
Check that the body and outlet of the instrument are clean and tidy, whether there are hairs, dust and fibers and other contaminants, if so, please carefully wipe clean.After you've been using the product for a while, be careful to see if your skin color has changed significantly —just in a bad direction, such as white and black. Purple, etc., you need to go to the hospital immediately for skin testing.
If you use the product for the first time, your skin color is within the available range: turn the gear in one gear. If there is no obvious pain, continue to increase the gear scale until the exposed skin has a slight pain. This is the most effective light gear.

Advantage of IPL laser hair remover

Remove excess hair and can permanently reduce hair growth, can effectively remove excess body hair, underarms. Bikini lines, arms and legs and other ideal areas of the body.

Post-care for IPL hair

Skin care

1. It is recommended to avoid touching water within 6 hours of lighting operation. The irradiated area can be moisturized with moisturizing water after 6 hours
2. Do not use perfume, masks and cosmetics containing matte beads after skin light
3. Please avoid UV exposure within 48 hours after lighting operation, and if necessary, plsapply sunscreen
4. Consult a surgeon if you have discomfort
5. Sunbathing is likely to cause the skin to become black or sunburn, at this time exposure will cause redness, inflammation and other conditions. Therefore, during the use of this product, pay attention not to sun exposure, sunburn repair time is not less than one week
6. In the light operation within 48 hours after the sun exposure area is strictly prohibited, even after 48 hours need to be exposed to the sun, but also always pay attention to whether the skin has a bad condition, please use sunscreen after 72 hours for sun protection
7. Artificial tanning skin is not allowed during lighting operation
8. Consult a doctor if you need plastic surgery after light (only in the area of irradiation)

Possible side effects

If shaved and lighted after the operation, for a period of time you feel dry and itchy skin, please rest assured. This reaction is harmless to the body and doesn't take long to disappear, and you can apply a mild, non-irritating moisturizer.

The safety of IPL laser hair remover

IPL laser hair removal is based on the principle of safety first and IPL pulse technology provides excellent safety and low energy levels.
Strong pulsed light technology enables long-term hair removal, resulting in a small portion of this energy being used on other optical hair remove de-vice, and low-energy IPL laser hair removal reduces possible injury or complications and contributes to overall safety.
Laser hair remover protect your eyes and skin! The laser hair remover has a built-in eye safety feature. It is designed to skin contact surface to touch the skin to flash work, laser hair removal device after 10 minutes do not operate, the device automatically shut down to stop working.

The safety of IPL laser hair remover

Please read important security information before you use it!
Laser hair remover are not designed for everyone. Please read and consider the values in the original section below before using it.
1. If you are not sure whether the laser hair remover is suitable for you, consult your surgeon or dermatologist.
2. Laser hair remover is a certain wavelength of strong pulsed light through the epidermis, into the dermis, selectively absorbed by the hair and hair follicles of the melanin particles, resulting in the thermal effect, the heat inside the hair conduction to the surrounding, hair follicles and stem cells, such as hair root complete destruction, resulting in permanent hair removal. Innate black skin tone contains higher melanin, pulsed bright light through the hair pigment can lead to burns on the black skin tone, bubbles and skin tone changes (pigmentation or pigmentation reduction), there is a high risk, this laser hair removal is not suitable for black skin color.

3. Do not use a laser hair remover on tanned skin or after recent exposure. Tanned skin, especially after exposure, contains a large amount of melanin, this skin type or color, containing those that do not seem to be tanned quickly, the presence of a large amount of melanin skin with laser hair removal device will have a higher risk. Use a laser
- 1.Power switch: power-on in 3 seconds, short press 1 second to switch gear.
- 2.Manual fash mode, press flash key 0.5 second into this mode(default in manual mode).
- 3.Auto fash mode, press flash key 3 seconds into auto mode, short press again to turn off this mode. epilator for recent expo-sure to the skin. Side effects such as burns, blisters, and changes in skin color, (excessive pigmentation or hypopigmentation).

Do not use under below circumstances

Specific conditions do not use laser hair remover for the time being, if you have the following circumstances please do not use this equipment;

1. If you are pregnant or lactating.
 2. If in the last 4 weeks you have been exposed or have been artifiially skinned.
 3. If you have a tattoo or permanent epilepsy in the area to be treated with. .
 4. If you have dark spots, such as freckles, birthmarks, moles in the area to be treated.
 5. If you have eczema, psoriasis, injury, exposed wounds or active infections, such as lip blisters in the treatment area, waiting until the wound heals before using this device.
 6. If you have a history of progenitor, a known photosensitivity (photosensitivity), or medication to make the skin more sensitive, including non-seedling anti-inflammatory agents (e. g. aspirin, ibuprofen, acetaminophen), tetracycline, sneezing, anti-vitamin A acid, sulfony uremia, sulfon-amide, tresel-amine, fluorine urine sneezing, changchunal, alkaloid, gray jaundice, fruit acid, all-anti vitamin A acid, acne special effects or local vitaminA.
 7. If you have a skin disease caused by diabetes, such as another systemic or metabolic disease.
 8. If you are currently using fruit acid, beta diacid, all anti vitamin A acid, acne special effects drug
- The above is the reference time, after all areas of treatment, it's better to repeat 2-3 times, which can avoid some areas are not irradiated.

- (iso-wit A acid), sib with vi-A acid or di-chloric acid, to avoid touching water within 6 hours of lighting operation. The irradiated area can then be moisturized with moisturizing water.
9. If you have taken acne special effects in the last 6 months (isotretinoin)
 10. If you have had Hrmong compound sorority in the last 3 months.
 11. If you have experienced an outbreak of herpes in the area you are treating, at least you consult your physician or receive preventive therapy before using this device.
 12. If you have epilepsy.
 13. If you have one-to-one active implants, such as starters, incontinence devices, islet pumps, etc.
 14. If you have related polymorphic diseases, such as tonic ulcer, polymorphic light outbreaks, solar urticaria, lupus erythematosus, etc.
 15. If you have a history of skin cancer or potential skin malignancies.
 16. If you do not confirm that the laser epilator is safe for use, consult your physician or dermatologist.
 17. If you have had radiotherapy or chemotherapy in the last 3 months.
 18. If you have some conditions a physician says you are not suitable for hair removal therapy.

How to use the laser hair remover safely.

Choose the intensity level carefully, as the intensity level increases, so does the risk of side effects. Always start using from the low intensity level, raise one level next time if no feel any discomforting. Avoid post use complications. Cover birthmarks and tattoos before use. Cover dark brown or black spots, such as freckles, birthmarks, moles, before use. Do not look directly at the device when light radiates from the window.
Do not use this device on nipples and reproductive organs (both male and female)
Do not use this device where you wish to grow hair later or other purposes other than hair removal. Do not use flash in open areas, during treatment, make sure surface fully contact with the skin.
If the skin contact surface is too hot, remove the device from skin.
10. Do not use flammable liquids such as alcohol (including perfumes, disinfectants, or other substances containing alcohol) or clean the skin with propylene when used.
11. Use of this equipment may cause temporary pigment changes.
12. Please keep this device away from children.
13. Please do not use this device on or for children.

Laser depilation treatment

usage plans and effects.
A normal complete hair growth cycle may take 18-24 months during this period, and multiple usage may be required. The effect of hair removal varies personal, depending on the color of the body's hair.
1-3 weeks
2-3 times/week
Hair began to fall naturally with a large number of hairs
4-5 weeks
2 times/week
New hair falls and inhabit the growth of hair, gain slippery skin
6-7 weeks
1 time/week
Less hair and gain more slippery skin
After 8 weeks
No more hair and keep the slippery skin, and you can use it if you want

How to use the IPL laser hair remover

1. Remove the laser hair remover and other accessories from the box.
2. Plug the power cord into the laser hair remover.
3. Plug the adapter into the socket.
4. First use a scraper to scrape the hairs on the surface of the skin clean. .
5. Aim at the skin, the flash mouth and skin to maintain vertical, press the flash key to ease hair.

Step by step for using

1. Scrape the hair of the hair-removal area with a scraper before start.
2. wear professional protective when operating, mirrors included in this product to protect the eyes.
3. Connect the power cord to the device power port, the power port has an LED on.
4. Press power switch for 3 seconds, gear light is on, the device can be used normally, there will be functional instructions and gear indicators.

5. After normal power-on, the function file is the "hair removal" function, if you need to switch function, while pressing the "On/Off function" key and"light emitting button", you can switch to "hair removal" , "Skin rejuvenation" , "Acne treatment" function.
6. After the normal power-up, is a gear, if you need to increase gear please short press the * switch machine / function" key, you can switch gear in turn. .
7. The product is vertically close to the skin, wait for the last gear light flashing after pressing the, light emission button, you can start hair removal (the last gear indicator flashes, indicating that the machine has detected the skin and into the working state), in the same place of the skin is best not more than 3 times to continuous flash hair removal.
8. Press and press the flash key 2 seconds to switch to automatic hair removal mode, automatic hair removal mode will automatically flash hair removal.

Reference time for depilation area

- 1.Power switch: power-on in 3 seconds, short press 1 second to switch gear.
- 2.Manual flash mode, press flash key 0.5 second into this mode(default in manual mode).
- 3.Auto flash mode, press flash key 3 seconds into auto mode, short press again to turn off this mode.

location	duration	Flash mode suggest
Armpit	4-5mins/side	Manual flash mode
Bikini line	8-10mins/side	Manual flash mode
Lip side	2-3mins/side	Manual flash mode
Thighs	25-30mins/side	Auto flash mode
Calf	25-30mins/side	Auto flash mode
Arm	25-30mins/side	Auto flash mode
Back	35-40mins	Auto flash mode

The above is the reference time, after all areas of treatment, it's better to repeat 2-3 times, which can avoid some areas are not irradiated.