



The hand-held
Diamond dermabrasion beauty
instrument

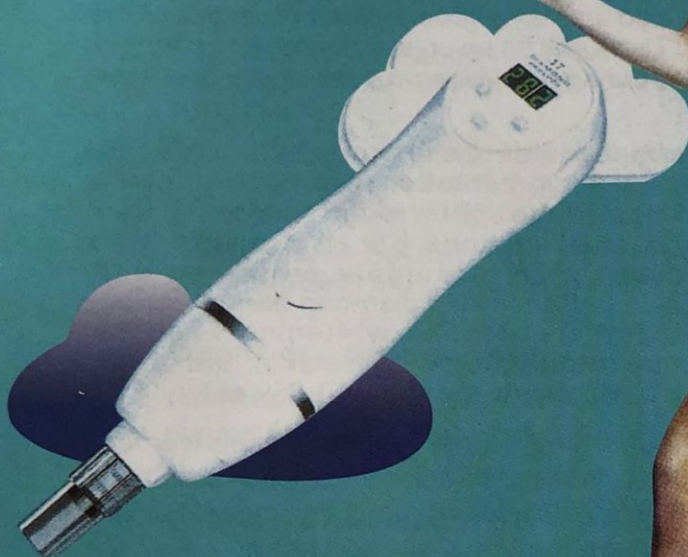

- Blackhead
- Diamond dermabrasion

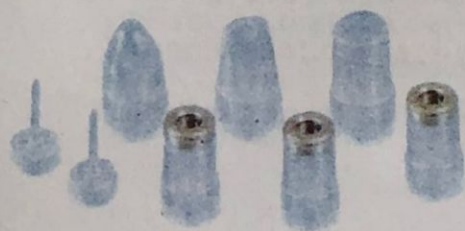

GCEC

Diamond dermabrasion beauty instrument introduction

diamond dermabrasion beauty instrument introduction

such as by the way. And with the vacuum suction apparatus, the suction microdiamond non-invasive, superficial dermabrasion physical nature.

This purely physical mechanical peel is based on the use of chemicals to remove skin, does not irritate skin with a mild advantage, appropriate. Combined for extremely sensitive skin. By suction, miniature size and thickness of the diamond. To control the depth of dermabrasion, less irritating, there will be no post-operation wound. Bothered to take care of, but also an easy pigmentation or redness of the skin problems, and there will be no crystal particles scatter the remaining questions during the whole process is very dry. Clean and safe.

Diamond dermabrasion beauty instrument can be used to improve the skin pigmentation disorder, compared with shallow pock youth, age spots, rough skin aging after sun exposure and facial fine lines, but also can promote skin renewal, immediately after surgery, presented with smooth skin. Touch... after making diamond dermabrasion, and then do the import or ultrasonic whitening. Since the United States to use formulation with ice skin, the effect will be more ideal, it will reach your beautiful wish!

In addition to diamond dermabrasion, this instrument is also equipped with a dedicated head of black smoke, ultra strong suction can be effectively sucked out pimples, blackheads and excess skin and fat, keep the skin clean and fresh.

Diamond dermabrasion function analysis

1. Damage

Microcirculation nutrition

Promote metabolism

Improve oxygen supply of nutrients

After the rapid removal of metabolic waste substances export

Lymphatic circulation (passive drainage system)

So that the cycle smoothly.

Eliminate excess water and old waste.

Increase nutrient supply function.

Adenoids (sebaceous glands, hair follicles)

So that the tube gland smooth clear obstruction.

Normal secretion of sebum nourish the skin.

● The treatment of acne, pock skin.

1. In-depth clear pores around the accumulation of sebum to prevent pore blockage, and strengthen the removal. Pock and surface dead cells.
2. To stimulate metabolism, promote oil and bacteria, toxins and old dead cells excluded, by. Plus the rate of new cell treatment pits;
3. Increased cell permeability and cell absorption, so that the active ingredient into the skin, inhibit fine. System of bacteria to heal the skin to appease the wound;
4. To promote epidermal cell rejuvenation, desalination pock, update skin, improve pock;

● On the dry, dry skin occurs.

1. The removal of aging dead skin cells, to prevent excessive accumulation of keratinocytes to keep the skin smooth, meticulous;
2. To promote blood, lymph circulation, so that water nutrients transfer more quickly, increase the brightness and luster;
3. Moisturizing repair effect: increase cell permeability, enhance product absorption, improve skin. Moisturizing and Moisturizing:
4. To promote epidermal basal cell regeneration, accelerate skin renewal, so that the skin tissue is restored. Often, let the transfer of water capacity to balance.

Treatment of relaxation, wrinkled skin occurs.

1. Accelerate the removal of the surface of the skin easy to damage the cell dust pollution, free radicals, enhance the cell's natural defense, to avoid collagen fibril degeneration.
2. To improve blood and lymph circulation, promote metabolism, so that the dermis to obtain adequate. Nutritional supplement, so that fine lines gradually disappear.

- 3、 Increase the permeability of the skin, so that more effective products into the deep, to promote new collagen fibers. Health, smooth wrinkles.
- 4、 To promote collagen fiber fibroblasts to produce collagen fibers, and increase skin defense function, To avoid free radicals and environmental damage, to keep skin compact and flexible.

● On dull, uneven color skin treatment.

- 1、 Remove the old dead skin cells, to avoid dead cell accumulation, keep the skin transparent and detailed.

2 100 content, so that skin rosy.

3 and moisture degree, so that the skin bright and shiny

4 To promote epidermal basal cell regeneration and accelerate skin renewal showing uniform, soft white new.

● The treatment of melanin calm skin.

- 1 To accelerate the removal of melanin in the epidermal keratinocytes and residual easily lead to the skin. Cell damage and environmental pollution and lipid peroxide substances can be added to accelerate the rapid dilution of embryonic black. Pigment spots.
- 2 To improve skin microcirculation, promote melanin metabolism, strengthen the deep beauty. White, reduction of melanin, the effect of faster and more obvious blemish.
- 3 Increase the permeability of the skin, so that whitening products in-depth cell absorption, effectively inhibit the black. The formation of elements to avoid melanin appear in the face.
- 4 To promote epidermal basal cell regeneration, and increase skin immunity strengthen the fight against UV lines, keep the skin white, bright.

Treatment effect

- 1 The hole 5-10 times the effect of 80%-100%
- 2 Large pores 1-2 times the effect of 90%
- 3 Atresia 2 times the effect of acne 100%
- 4 Lposuction dermatitis 3-5 times the effect of 100%
- 5 Tough skin 1-2 times the effect of 100%
- 6 Fiber lump 2-3 times the effect of 100%

7 nasolabial wrinkles effect as the case of 80%-100%

8 the projection effect of age spots, as the case of 80%-100%

9 the effect of small wrinkles and rosacea, as the case of 80%-100%

10 tight lift 1-10 times the effect of 80%-100%

11 pock scars and 3-5 times the effect of 90%

Handheld Diamond Dermabrasion Beauty Instruments

● Small, lightweight, easy to carry;

Host only 160g, home, easy to carry out.

● 3 different degrees of thickness miniature head;

Fine carved diamond grit drill head for relatively delicate and sensitive parts of the skin (such as eye periphery): thickness diamond grit drill ships carved heads for the face, neck: drilling Coarse stone carving drill heads for rough skin and wrinkles more parts (eg body, elbows, knees).

● 3 different sizes of black smoke treatment head; 1

Fine diamond grit drill carved heads for delicate skin or relatively sensitive parts (such as around the eyes).

● Double filtration device, easy to clean, and effectively prevent the blockage into the inside of the machine.

● For a variety of skin, sensitive skin can also be done.

● Non-irritating, non-stinging, no need anesthesia.


● In physical mechanical dermabrasion can be easily controlled, no side effects, safe and reliable.

● Non-invasive treatment it is not irritating, no wound.

● You can see the immediate effect immediately.

● After completion of care will not affect the work and socializing.

Diamonddermabrasion beaulynstrument
Installation instructions

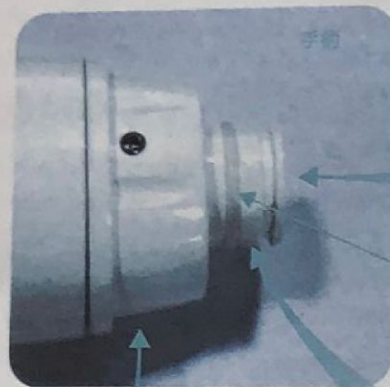


the first step:
Installation sponge filter

sponge



Silicone ring

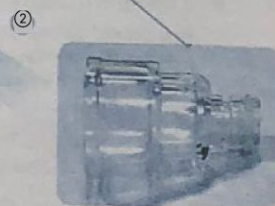

the Step two:
Install a transparentconnection




the third step:
Install aprotectivefrontofthehandlehousing



Rotarymodeinstallation,completeinspection
Inner siliconering is still in the slot,ifhair
Now is extruded silicone ring slot need to
reinstall Loading.




the fourth step:
Installation of diamond derma
brasion head or black smoke
Treatment head

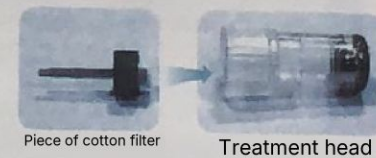

Piece of cotton filter

Treatment head

the fifth step:
Install a transparent
protective sleeve

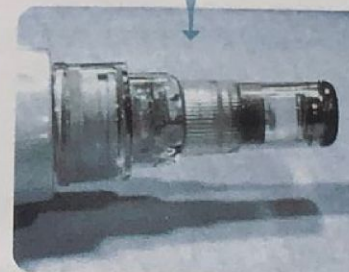


Handle with DC12V power supply.
When not using the available cover
dust cover The interface cover.

Installation diagram



- Intact siliconing various locations to ensure the tightness of the entire system. Do not pay attention when installing the missing equipment, check the silicone ring is in good condition when installed.
- Use 3-5 times the need to replace a piece of cotton piece of cotton inside the filter
- Regular inspection sponge filter is blocked, if found a large obstruction. Please change the sponge material

Operation Panel

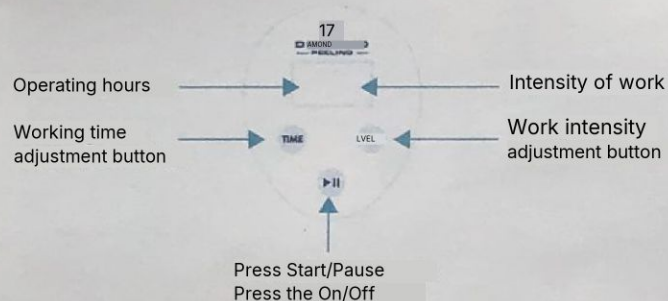

Instructions

Face operating practices




- 1, Cleaning the face, wipe clean with a wet towel, the face cannot be too dry, so as not to affect the results,
- 2, Long press the start/pause button to open the instrument, adjust the working time and work intensity;
- 3, The handle treatment head close to the skin so that it gives pressure to the skin. Inhalation treatment head, hold hands. The handle moves back and forth on the skin;

Precautions after use

- 1, Regular use of office or live cells or collagen products used for repair and stabilization.
- 2, Be sure to do sun protection work.
- 3, For those who apply makeup, please use pure mineral powder to avoid exclusion.
- 4, Please wash your face with cold water.
- 5, 10-15 days after treatment do not wash with warm water or contact heat.
- 6, Due to personal skin condition, vitality and treatment of different shades of about 1-3 weeks interval. In order to carry out treatment.

Bill of Material

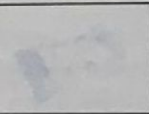Serial number	Pictures	Name
1		Machine
2		Power supply
3		Cotton bag filter.
4		Filtration with sponge rubber.
5		Bag
6		Manual
7		Masonry slightly carved.
8		Blackhead acne.
9		Filter assembly.

Service and Repair

1. This machine white quality problems from the date of purchase (with tickets) within a week and look good before a new machine, the host Paul Repair for one year. other onwaatoutofwaatyiodalsoesoileoteman and updating, recovery fee, spare parts And consumer goods.

2. ellou

- ① Disassemble, modify the product caused by the failure;
- ② Careless use of the beat, drop caused by discharge barrier;
- ③ Due to lack of reasonable care caused by failure;

Maintenance Record Form

Product Name:		Purchase Date:	
Dealership Name:		Phone:	
Address:			
Customer Name:		Phone:	
Address:			
Maintenance Date	Maintenance Content	cause of issue	