



The hand-held
Diamond dermabrasion beauty
instrument

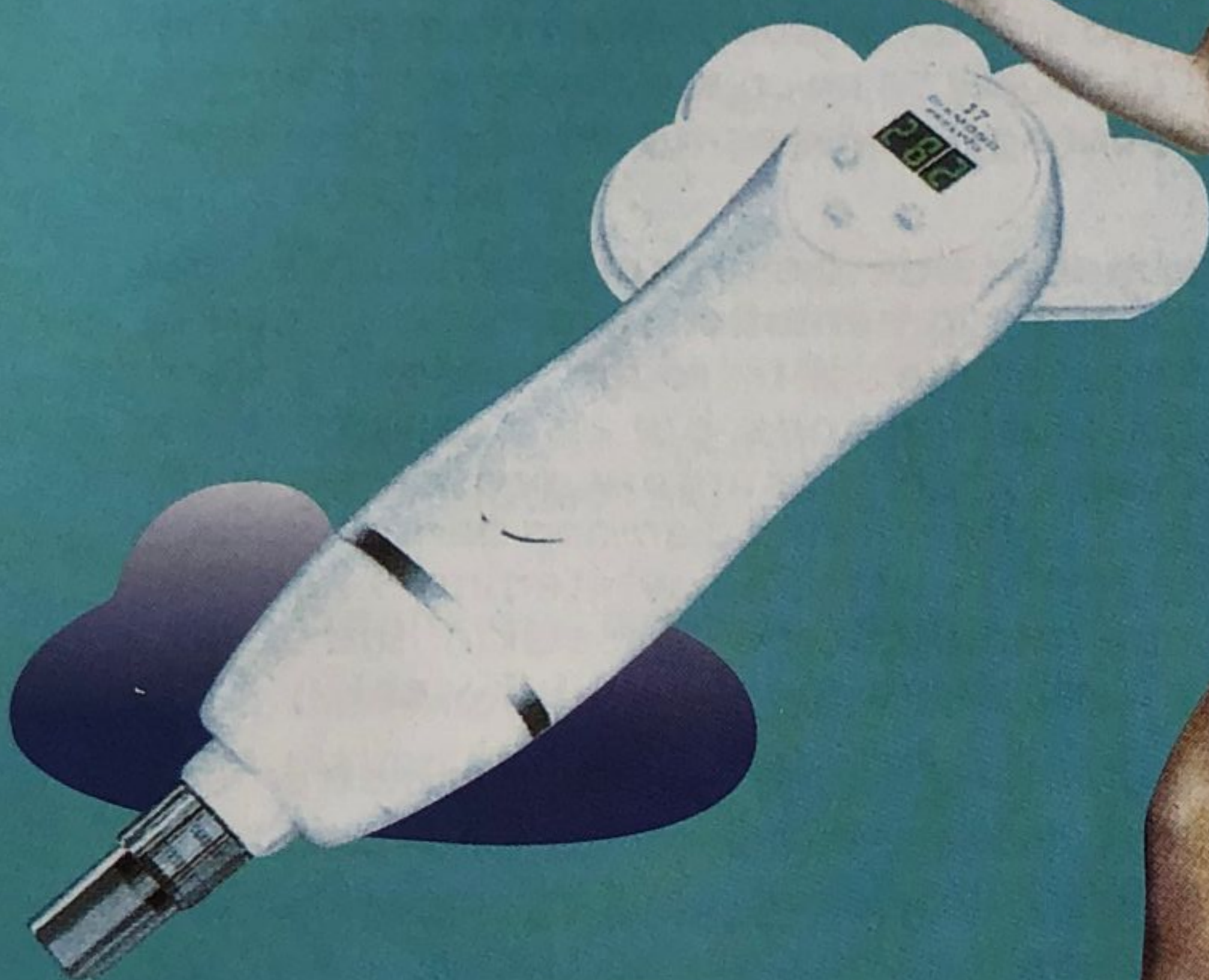

- Blackhead
- Diamond dermabrasion

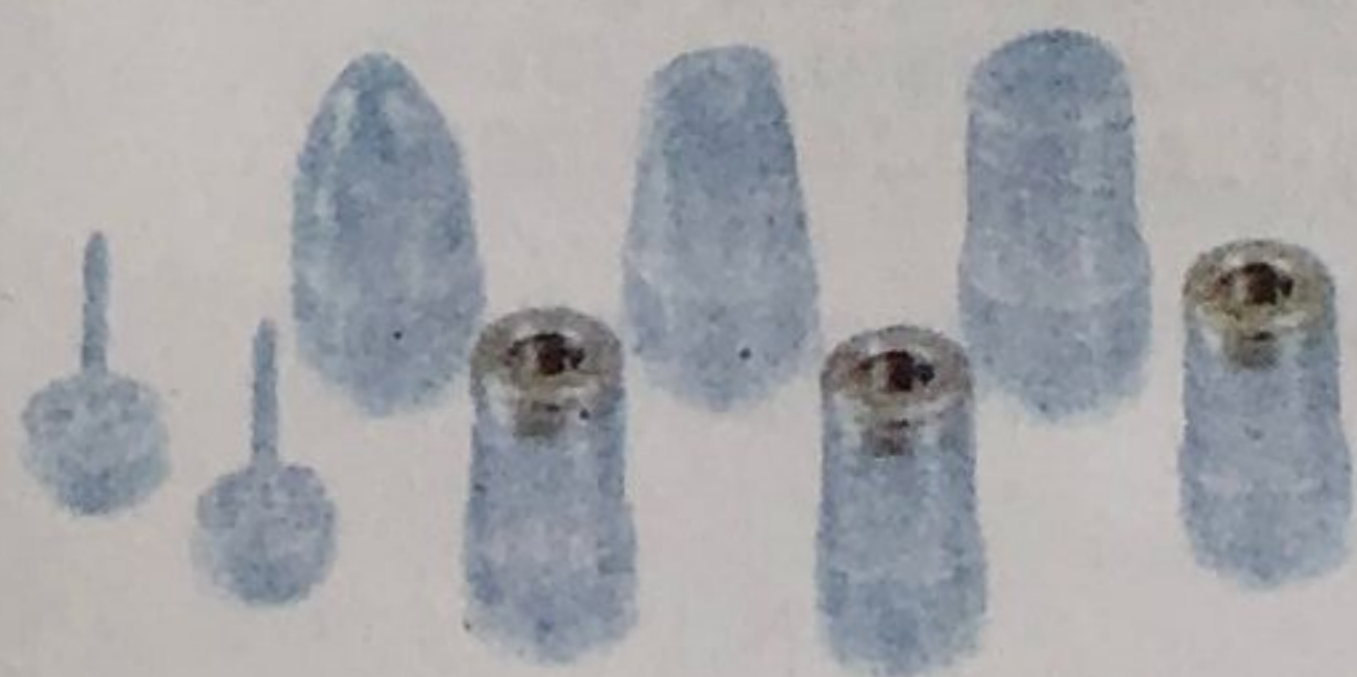

Diamond dermabrasion beauty instrument Introduction

Diamond dermabrasion beauty is the use of patented diamond dermabrasion tube, combined with inlaid different thickness of the diamond grit drill carved head, rubbing back and forth by the way. And with the vacuum suction apparatus suction control the intensity of the skin shallow grind. Off layer is a non-invasive, superficial dermabrasion physical nature.

This purely physical, mechanical principles, compared with acid peels be based on Way chemicals to remove skin, does not irritate skin with a mild advantage, appropriate Combined for extremely sensitive skin. By suction miniature size and thickness of the diamond To control the depth of dermabrasion, less irritating, there will be no post-operation wound Bothered to take care of, but also an easy pigmentation or redness of the skin problems, and There will be no crystal particles scatter the remaining questions during the whole process is very dry Clean and safe.

Diamond dermabrasion beauty instrument can be used to improve the skin pigmentation disorder, compared with Shallow pock youth, age spots, rough skin aging after sun exposure and facial Fine lines, but also can promote skin renewal, immediately after surgery, presented with smooth skin Touch after making diamond dermabrasion, and then do the import or ultrasonic whitening Since the United States to use formulation with ice skin, the effect will be more ideal, it will reach Your beautiful wish!

In addition to diamond dermabrasion, this instrument is also equipped with a dedicated head of black smoke, ultra Strong suction can be effectively sucked out pimples, blackheads and excess skin and fat, Keep the skin clean and fresh.

Diamond dermabrasion function analysis

- Diamond dermabrasion role in the activation of the dermis to make, repair, regeneration

Microcirculation (nutrition, waste exchange places)

Promote metabolism

Improve oxygen supply of nutrients

After the rapid removal of metabolic aging substance export

Lymphatic circulation (passive drainage system)

So that the cycle smoothly.

Eliminate excess water and old waste.

Increase nutrient supply function.

Adrenal glands (sebaceous glands, hair follicles)

So that the tube gland smooth clear obstruction.

Normal secretion of sebum nourish the skin.

● The treatment of acne, pock skin.

- 1、In-depth clear pores around the accumulation of sebum to prevent pore blockage, and strengthen the removal. Pock and surface dead cells.
- 2、to stimulate metabolism, promote oil and bacteria, toxins and old dead cells excluded, by. Plus the rate of new cell treatment pits;
- 3、Increased cell permeability and cell absorption, so that the active ingredient into the skin, inhibit fine. System of bacteria to heal the skin to appease the wound;
- 4、to promote epidermal cell rejuvenation, desalination pock, update skin, improve pock;

● On the dry, dry skin occurs.

- 1、the removal of aging dead skin cells, to prevent excessive accumulation of keratinocytes to keep the skin smooth, meticulous;
- 2、promote blood, lymph circulation, so that water nutrients transfer more quickly, increase the brightness and. luster;
- 3、moisturizing repair effect: increase cell permeability, enhance product absorption, improve skin. Moisturizing and Moisturizing:
- 4、to promote epidermal basal cell regeneration, accelerate skin renewal, so that the skin tissue is restored. Often, let the transfer of water capacity to balance.

Treatment of relaxation, wrinkled skin occurs.

- 1、accelerate the removal of the surface of the skin easy to damage the cell dust pollution, free radicals, Enhance the cell's natural defense, to avoid collagen fibril degeneration.
- 2、To improve blood and lymph circulation, promote metabolism, so that the dermis to obtain adequate. Nutritional supplement, so that fine lines gradually disappear.

- 3、 Increase the permeability of the skin, so that more effective products into the deep, to promote new collagen fibers. Health, smooth wrinkles.
- 4、 To promote collagen fiber fibroblasts to produce collagen fibers, and increase skin defense function, To avoid free radicals and environmental damage, to keep skin compact and flexible.

● On dull, uneven color skin treatment.

- 1、 Remove the old dead skin cells, to avoid dead cell accumulation, keep the skin transparent and detailed.
- 2、 Improve blood, lymph circulation, promote metabolism, increase skin oxygen content, so that. Skin rosy.
- 3、 Increase cell permeability, enhance product absorption, improve skin insulation and moisture degree, So that the skin bright and shiny.
- 4、 To promote epidermal basal cell regeneration, accelerate skin renewal, showing uniform, soft white new.

● The treatment of melanin calm skin.

- 1、 To accelerate the removal of melanin in the epidermal keratinocytes and residual easily lead to the skin. Cell damage and environmental pollution and lipid peroxide substances can be added to accelerate the rapid dilution of embryos black. Pigment spots.
- 2、 To improve skin microcirculation, promote melanin metabolism, strengthen the deep beauty. White, reduction of melanin, the effect of faster and more obvious Blemish.
- 3、 Increase the permeability of the skin, so that whitening products in-depth cell absorption, effectively inhibit the black. The formation of elements to avoid melanin appear in the face.
- 4、 To promote epidermal basal cell regeneration, and increase skin immunity, strengthen the fight against UV. Line, keep the skin white, bright.

Treatment effect


- 1、 The hole 5 - 10 times the effect of 80% - 100%
- 2、 Large pores 1-2 times the effect of 90%
- 3、 Atresia 2 times the effect of acne 100%
- 4、 Lposuction dermatitis 3-5 times the effect of 100%
- 5、 Tough skin 1-2 times the effect of 100%
- 6、 Fiber lump 2-3 times the effect of 100%

- 7、 nasolabial wrinkles effect as the case of 80% - 100%
- 8、 the projection effect of age spots, as the case of 80% - 100%
- 9、 the effect of small wrinkles and rosacea, as the case of 80% - 100%
- 10、 tight lift 1-10 times the effect of 80% - 100%
- 11、 pock scars and 3-5 times the effect of 90%

The hand-held diamond dermabrasion beauty instrument features

- Small, lightweight, easy to carry;
Host only 160g, home, easy to carry out.
- 3 different degrees of thickness miniature head;
Fine carved diamond grit drill head for relatively delicate and sensitive parts of the skin (such as eye Eye peripheral) ; thickness diamond grit drill ships carved heads for the face, neck: drilling Coarse stone carving drill heads for rough skin and wrinkles more parts (eg body Body, elbows, knees).
- 3 different sizes of black smoke treatment head;
Fine diamond grit drill carved heads for delicate skin or relatively sensitive parts (such as Around the eyes).
- Double filtration device, easy to clean, and effectively prevent the blockage into the inside of the machine.
- For a variety of skin, sensitive skin can also be done.
- Non-irritating, non-stinging, no need anesthesia.
- In physical mechanical dermabrasion can be easily controlled, no side effects, Safe and reliable.
- Non-invasive treatment, it is not irritating, no wound.
- You can see the immediate effect immediately.
- After completion of care will not affect the work and socializing.

Diamond dermabrasion beauty instrument installation instructions

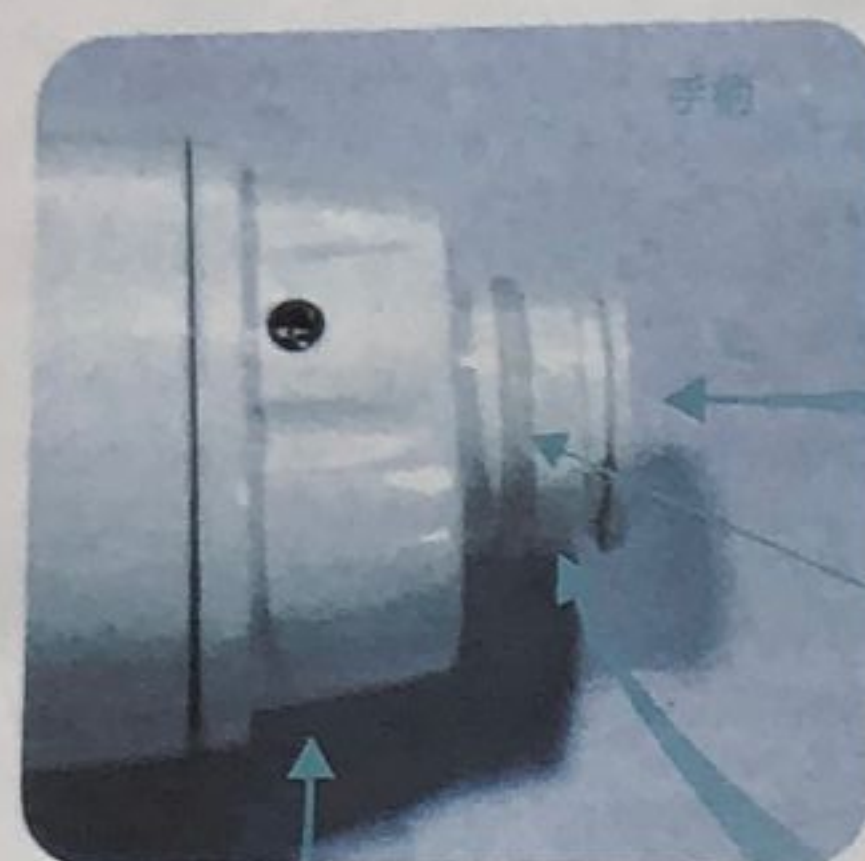

the first step:
Installation sponge filter

sponge

①



Sponge filter

Silicone ring

②



the Step two:
Install a transparent connection



the third step:
Install a protective front of the handle housing



Rotary mode installation, complete inspection
Inner silicone ring is still in the slot, if hair
Now is extruded silicone ring slot need to
reinstall Loading.



the fourth step:
Installation diamond derma
brasion head or black smoke
Treatment head

④




Piece of cotton filter



Treatment head

the fifth step:
Install a transparent
protective sleeve

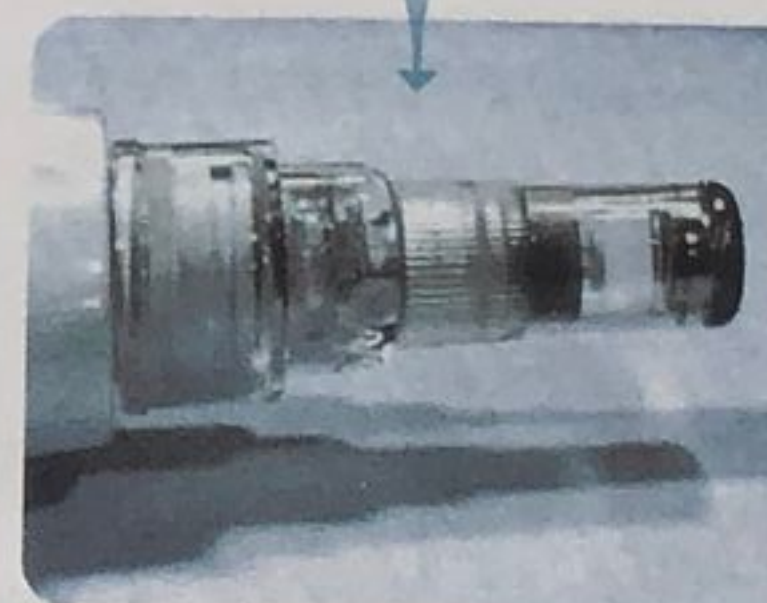


⑤




Handle with DC12V power supply,
When not using the available cover
dust cover The interface cover.

Installation diagram

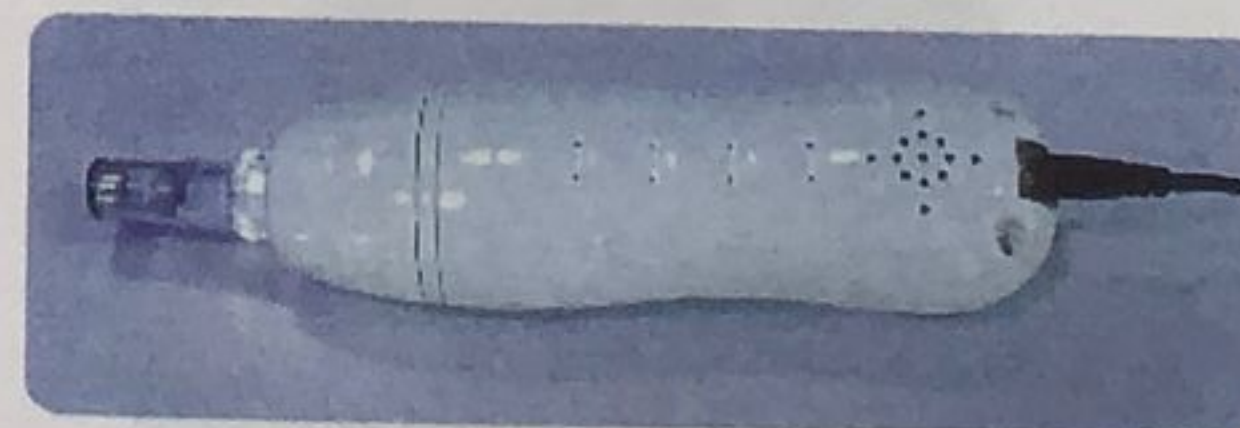

- Intact silicone ring various locations to ensure the tightness of the entire system, Do not pay attention when installing the missing equipment, check the silicone ring is in good condition when installed.
- Use 3-5 times the need to replace a piece of cotton piece of cotton inside the filter
- Regular inspection sponge filter is blocked, if found a large obstruction Please change the sponge material

Operation Panel

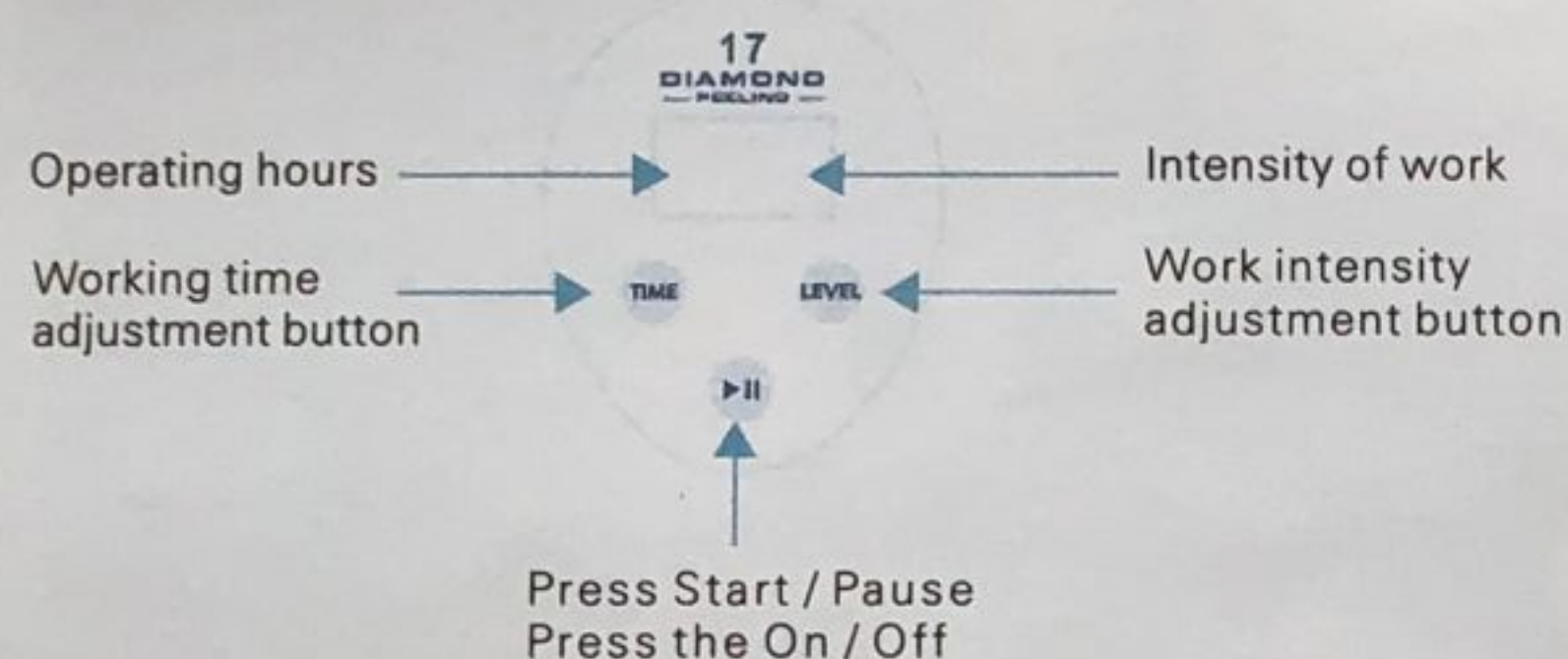

Instructions

Face operating practices





- 1, Cleaning the face, wipe clean with a wet towel, the face can not be too dry, so as not to affect the results,
- 2, Long press the start / pause button to open the instrument, adjust the working time and work intensity;
- 3, The handle treatment head close to the skin, so that negative pressure to the skin inhalation treatment head, hold hands. The handle moves back and forth on the skin;

Precautions after use

- 1、 Regular use of ice or live cells or collagen products used for repair and stabilization.
- 2、 Be sure to do sun protection work.
- 3、 For those who apply makeup, please use pure mineral powder to avoid exclusion.
- 4、 Please wash your face with cold water.
- 5、 10-15 days after treatment do not wash three warm or contact heat.
- 6、 Due to personal skin condition, vitality and treatment of different shades of about 1-3 weeks interval. In order to carry out treatment.

Bill of Material

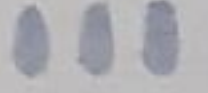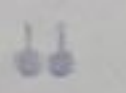Serial number	Pictures	Name
1		Machine
2		Power supply
3		Cotton bag filter.
4		Filtration with sponge rubber.
5		Bag
6		Manual
7		Masonry slightly carved.
8		Blackhead acne.
9		Filter assembly.

Service and Repair

1. This machine white quality problems from the date of purchase (with tickets) within a week and look good before a new machine, the host Paul Repair for one year, other non - warranty or out of warranty period, also responsible for the maintenance and updating, recovery fee, spare parts And consumer goods.

2. The following users' personal reasons for failure, not provide free warranty service:

- ① Disassemble, modify the product caused by the failure;
- ② Careless use of the beat, drop caused by discharge barrier;
- ③ Due to lack of reasonable care caused by failure;

Maintenance Record Form

Product Name:

Purchase Date:

Dealership Name:

Phone:

Address:

Customer Name:

Phone:

Address:

Maintenance Date

Maintenance Content

cause of issue